

Ventaire™ Scope Drying and Tracking Retrofit Unit




Support scope drying and storage processes, safety protocols, and reduce the chances (and time and expense) of reprocessing expired endoscopes with the Ventaire scope drying and tracking system. The system meets the 2021 Multisociety Task Force Scope Storage and Drying Guidelines.

The Ventaire system continuously filters air into internal scope channels to efficiently and thoroughly dry scopes. The easy-to-operate system features timed drying and venting cycles. HEPA filtered air removes 99.9% retention rate of particles less than 0.1µm. It works with all major endoscope brands.

- Continuous venting of HEPA filtered air through endoscope channels, providing safe scope storage 24/7
- 10-minute default drying cycle
- Modify or add to scope inventory at any time
- Customizable report captures all scope information and status
- Monitor scope storage from touchscreen on cabinet or remotely via website (subscription required)

Ventaire Scope Drying and Tracking Retrofit Unit

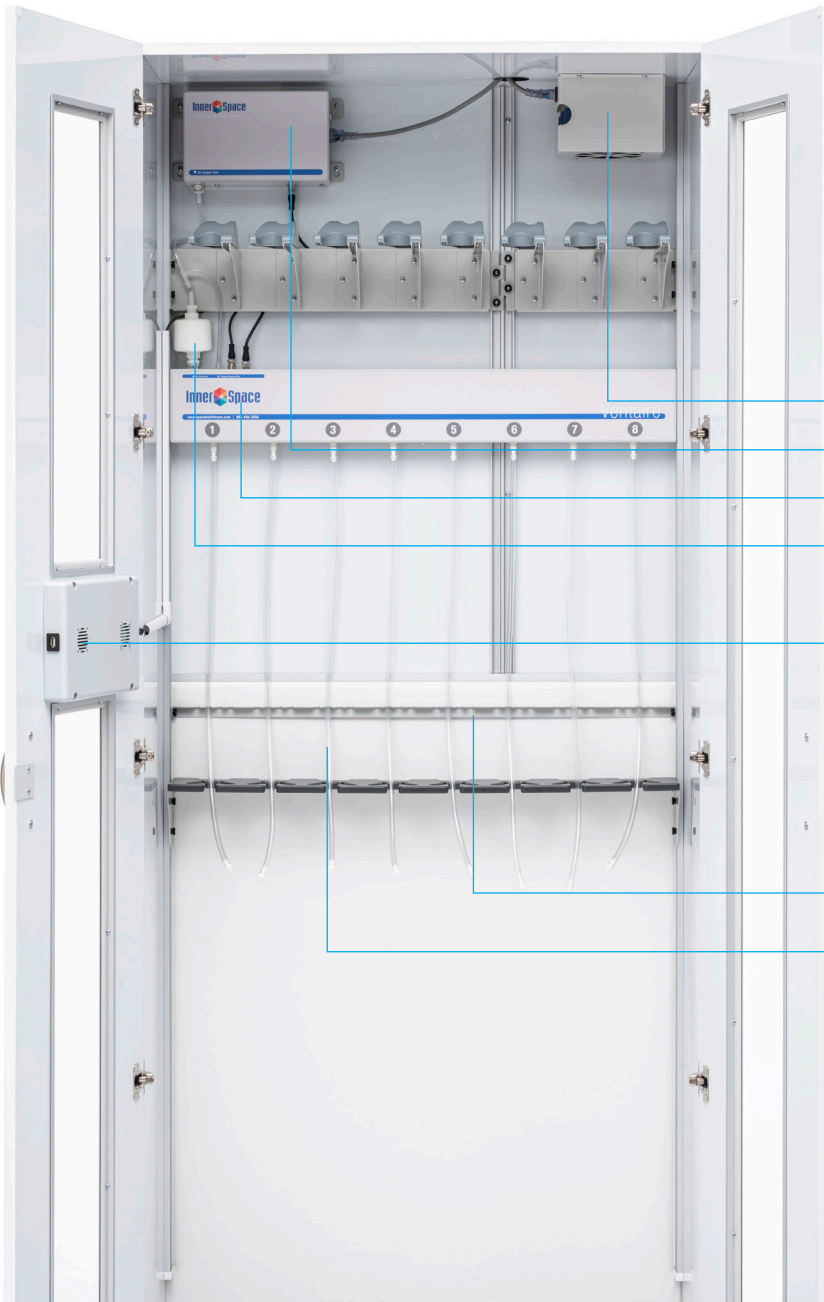

The retrofit unit installs in existing Evolve scope cabinets with adjustable scope managers and cord retainers. It includes a Ventaire scope drying and tracking unit, filters and tubes, and replacement door with 10" touchscreen. The unit is designed for easy installation and start-up.

Optional HEPA fan unit provides positive-pressure air in cabinet

Energy-efficient, quiet air compressor flows air into manifold

Manifold calibrates air circulation into active scope channels

Filter controls air circulation into endoscopes

Easy-to-operate 10" touchscreen



Tube management keeps endoscope tubes organized and secure

Tubes connect to components in the Ventaire system to vent HEPA filtered air through endoscope channels

Ventaire Scope Drying and Tracking Retrofit Unit

Dimensions

Compressor: 9.61"w x 3.56"d x 5.96"h

5- and 10-port manifold: 24.30"w x 0.64"d x 4.88"h

8- and 16-port manifold: 33.62"w x 0.64"d x 4.88"h

9- and 18-port manifold: 37.6"w x 0.64"d x 4.88"h

Interior Configuration:

26.75"w cabinet: holds a 5- or 10-port manifold

36"w cabinet: holds a 8- or 16-port manifold

40"w cabinet: holds a 9- or 18-port manifold

Specifications

Compressor air output:

- up to 45 LPM

Compressor decibels:

- 62 dB, with doors closed


Electrical requirements for tracking unit:

- 120 VAC, 60Hz, 15-amp standard duplex receptacle
- Includes a 125V, 10-amp, 10-foot-long power cord

Options

- Subscription to InnerSpace Smart Solutions website for remote monitoring of scopes

Product Number Reference

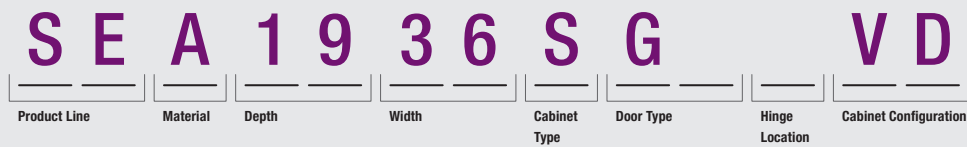

Product Numbers

To order a Ventaire scope drying and tracking retrofit unit, first specify the replacement door size to match the cabinet size, and then specify the drying and tracking retrofit unit.

Door for 26.75" w x 48" h left-hinged cabinet

- AireCore: SEA1927SGLWVD (glass door), SEA1927SSLWVD (solid door)
- AireCore with Brushed Aluminum: SEB1927SGLWVD (glass door), SEB1927SSLWVD (solid door)
- Stainless steel: SES1927SGLWVD (glass door), SES1927SSLWVD (solid door)

Door for 26.75" w x 48" h right-hinged cabinet

- AireCore: SEA1927SGRWVD (glass door), SEA1927SSRWVD (solid door)
- AireCore with Brushed Aluminum: SEB1927SGRWVD (glass door), SEB1927SSRWVD (solid door)
- Stainless steel: SES1927SGRWVD (glass door), SES1927SSRWVD (solid door)

Door for 26.75" w x 60" h left-hinged cabinet

- AireCore: SEA1927SGLBVD (glass door), SEA1927SSLBVD (solid door)
- AireCore with Brushed Aluminum: SEB1927SGLBVD (glass door), SEB1927SSLBVD (solid door)
- Stainless steel: SES1927SGLBVD (glass door), SES1927SSLBVD (solid door)

Door for 26.75" w x 60" h right-hinged cabinet

- AireCore: SEA1927SGRBVD (glass door), SEA1927SSRBVD (solid door)
- AireCore with Brushed Aluminum: SEB1927SGRBVD (glass door), SEB1927SSRBVD (solid door)
- Stainless steel: SES1927SGRBVD (glass door), SES1927SSRBVD (solid door)


Door for 26.75" w x 92" h left-hinged cabinet

- AireCore: SEA1927SGLVD (glass door), SEA1927SSLVD (solid door)
- AireCore with Brushed Aluminum: SEB1927SGLVD (glass door), SEB1927SSLVD (solid door)
- Stainless steel: SES1927SGLVD (glass door), SES1927SSLVD (solid door)

Door for 26.75" w x 92" h right-hinged cabinet

- AireCore: SEA1927SGRVD (glass door); SEA1927SSRVD (solid door)
- AireCore with Brushed Aluminum: SEB1927SGRVD (glass door); SEB1927SSRVD (solid door)
- Stainless steel: SES1927SGRVD (glass door); SES1927SSRVD (solid door)

Product Number Reference

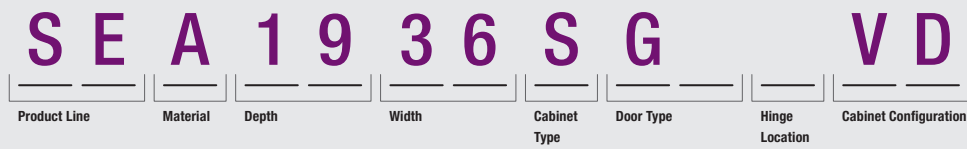

Product Numbers

To order a Ventaire scope drying and tracking retrofit unit, first specify the replacement door size to match the cabinet size, and then specify the drying and tracking retrofit unit.

Left-hinged door for 36" w cabinet

AireCore: SEA1936SGLVD (glass door); SEA1936SSLVD (solid door)

AireCore with Brushed Aluminum: SEB1936SGLVD (glass door); SEB1936SSLVD (solid door)

Stainless steel: SES1936SGLVD (glass door); SES1936SSLVD (solid door)

Left-hinged door for 40" w cabinet

AireCore: SEA1940SGLVD (glass door); SEA1940SSLVD (solid door)

AireCore with Brushed Aluminum: SEB1940SGLVD (glass door); SEB1940SSLVD (solid door)

Stainless steel: SES1940SGLVD (glass door); SES1940SSLVD (solid door)

Product Number Reference



Product Numbers

Ventaire Scope Drying and Tracking Retrofit Unit – For 26.75" w left-hinged cabinet

5-port unit: SVDURL5

10-port unit: SVDURL10

Ventaire Scope Drying and Tracking Retrofit Unit – For 26.75" w right-hinged cabinet

5-port unit: SVDURR5

10-port unit: SVDURR10

Ventaire Scope Drying and Tracking Retrofit Unit – For 36" w cabinet

8-port unit: SVDUR8

16-port unit: SVDUR16

Ventaire Scope Drying and Tracking Retrofit Unit – For 40" w cabinet

9-port unit: SVDUR9

18-port unit: SVDUR18

Ventaire Scope Drying and Tracking Retrofit Unit

Replace these items every 6 months for optimal performance:

- Compressor tube/air filter replacement, package of 1: SVCTAF
- Manifold tube replacement, package of 5: SVMT5
- Manifold tube replacement, package of 8: SVMT8
- Manifold tube replacement, package of 9: SVMT9
- Manifold auxiliary water tube replacement, package of 5: SWWAT5
- Manifold auxiliary water tube replacement, package of 8: SWWAT8
- Manifold auxiliary water tube replacement, package of 9: SWWAT9
- HEPA filter replacement, package of 1: SHEPFIL, package of 3: SHEPFIL3

Single-use connector kits:

To find which connector kit is compatible with your endoscope, visit <https://www.innerspacehealthcare.com/Ventaire>

- Connector Kit 1, package of 100: SVCK1
Includes connector and biopsy port cover
- Connector Kit 2, package of 100: SVCK2
Includes connector, biopsy port cover, and auxiliary water connector
- Connector Kit 3, package of 100: SVCK3
Includes connector and luer plug
- Connector Kit 4, package of 100: SVCK4
Includes connector, water jet adaptor, and luer plug
- Connector Kit 5, package of 100: SVCK5
Includes connector, biopsy port cover, and luer coupler

To order, call an InnerSpace Customer Service representative at 888-435-2256 or visit <https://www.innerspacehealthcare.com/Ventaire>

Ventaire Scope Drying and Tracking Retrofit Unit

Materials and Finishes

Compressor:

Powder-coated steel

Manifold:

Powder-coated steel

Tubing:

Soft PVC

Compressor tube/air filter:

Soft PVC tube and polypropylene filter

Connectors (vary by type):

HDPE, nylon (white), polypropylene (white) and polycarbonate (clear)